




## REF. 1150 – CALLING CARD WET FILM GAUGE

Format complies with ASTM D4414 – Procedure A, BS3900 Part C5, ISO 2808 Method 7B

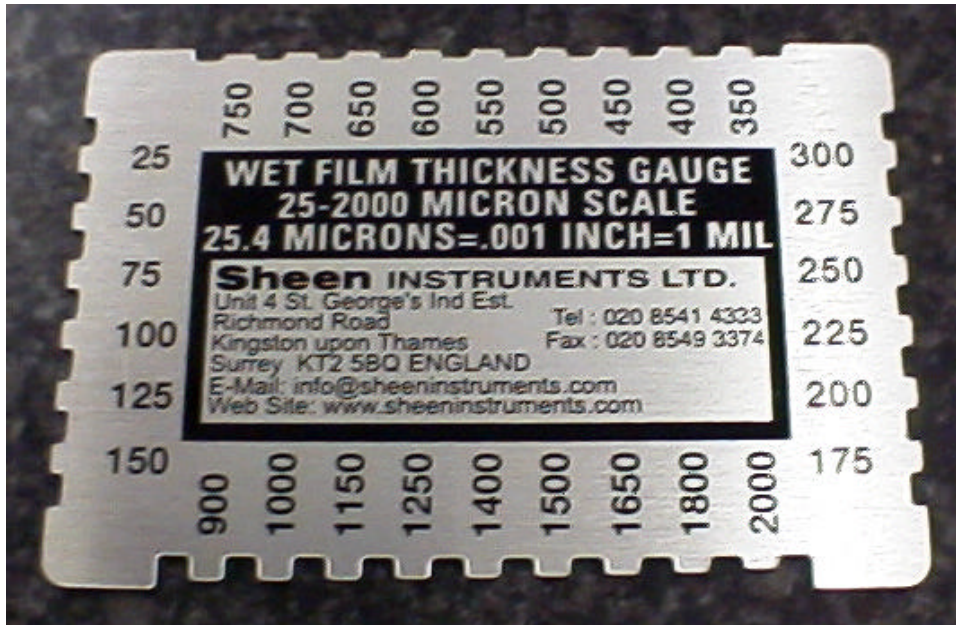

A precision made wet film thickness gauge to check wet film thicknesses from 25 – 2000 microns

### DIRECTIONS FOR USE

1. Select one of the four sides of the gauge according to the approximate wet film thickness, if known.
2. Hold the gauge perpendicular to surface and press it into the wet coating until the gauge rests squarely on the metal base ( substrate ).

DO NOT DRAG OR TWIST THE GAUGE

3. Withdraw the gauge and note the **deepest** tooth which has picked up the coating on it.

THE WET FILM THICKNESS LIES BETWEEN THE TWO READINGS.

4. After use clean immediately in appropriate solvent.

### IDEAL PROMOTIONAL ITEM

We can print your own company logo or name and address on the front or reverse side of the gauge - minimum quantity 100 gauges.

Particularly suitable for chemical and paint manufacturers, inspectors and useful for anyone who needs to know the wet film thickness of any coating material, adhesive, pipe or roof coating etc.

Size: 82 x 57mm ( as illustrated )  
Supplied in a transparent plastic folder.

Owing to continuous development, we reserve the right to introduce improvements and modify specifications without prior notice.

Sheen Instruments Ltd

Unit 4, St. Georges Ind. Est., Richmond Road, Kingston, KT2 5BQ England.

Tel: 020 8541 4333 Fax: 020 8549 3374 Intl.Tel: 44 20 8541 4333 Intl.Fax: 44 20 8549 3374

Internet : [www.sheeninstruments.com](http://www.sheeninstruments.com)

Email: [info@sheeninstruments.com](mailto:info@sheeninstruments.com)